


Acetabular Reamers

A smooth, quick ream for accurate and consistent sizing.

Designed and manufactured for performance, Enztec Acetabular Reamers provide a smooth, quick ream. Their ultra-sharp teeth, achieved through the use of an innovative photochemical machining process, result in accurate and consistent sizing. This allows the Enztec Acetabular Reamers to be used in single-shot reaming techniques and with lower cutting forces.

With a variety of sizes available, and Crossbar or Bridgeback connections, there is an Enztec Acetabular Reamer to suit your needs.

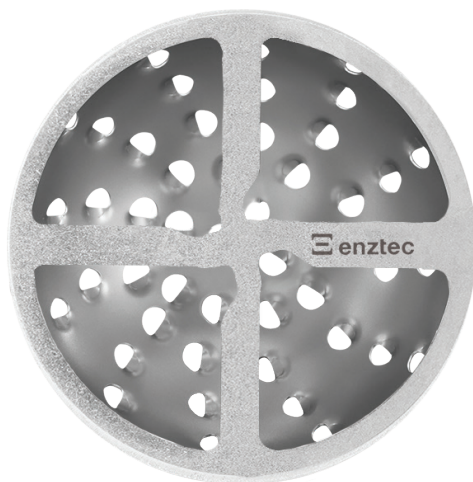

— CROSSBAR



— BRIDGEBACK

Features

- Sharp cutting edge achieved through photochemical machining process
- Accurate and consistent sizing
- Variety of sizes to support different patient anatomies (sizes range from 42 - 66)
- Crossbar or Bridgeback connections

Specifications

- Materials: Stainless Steel
- Weight: 22-42g / 0.05-0.10lbs (each)
- Diameter: 42mm to 66mm
- Sizing increments: 1mm